
VALO[®]

GRAND

CORDED BROADBAND **LED** CURING LIGHT

UP ULTRADENT
PRODUCTS, INC.

WARNINGS AND PRECAUTIONS

Read all instructions before operating this unit. The manufacturer accepts no liability for any damage resulting from the improper use of this unit and/or for any purpose other than those covered by these instructions.

- VALO GRAND Corded curing light emits high intensity light very similar to high intensity Quartz halogen lights and must only be used as indicated in this manual.
- DO NOT look directly into the light output. Patient, clinician, and assistants should always wear amber UV eye protection when VALO GRAND Corded curing light is in use.
- USE CAUTION when treating patients who suffer from adverse photobiological reactions or sensitivities, patients who are undergoing chemotherapy treatment, or patients being treated with photosensitizing medication.
- DO NOT expose soft oral tissues at close proximity for more than 10 seconds in any mode. If longer curing time is required, use multiple shorter curing cycles to avoid heating soft tissue, or use a dual-cure product.
- If using VALO GRAND Corded curing light in close proximity to the gingiva, DO NOT expose tissue for more than 20 seconds. If a 40-second cure is needed, allow 10 seconds between two 20-second cures. If longer curing is required, consider a dual-cure product.
- The Xtra Power mode has a 2-second safety delay to limit heating during consecutive curing. At the end of the delay, beeping indicates unit is ready for continued use.
- Take extra care to avoid directing light toward soft tissues.
- WARNING: No modification of this equipment is allowed.
- When disposing of electronic waste, (i.e. devices, chargers, batteries and power supplies), follow local waste and recycling guidelines.

General Information:

Indications for Use: The source of illumination for curing photo-activated dental restorative materials and adhesives.


VALO GRAND Corded curing light is an LED curing light for the polymerization of light cured dental materials. With its broadband spectrum, VALO GRAND Corded curing light is designed to polymerize all light cured products in the wavelength range of 385-515nm per ISO 10650.

VALO GRAND Corded curing light has a medical grade, international power supply and is suitable for power outlets from 100 to 240 volts. The handpiece is designed to rest in a standard dental unit bracket or can be custom mounted using the bracket included with the kit.

Product Components:

- 1 – VALO GRAND Corded curing light with 7 foot / 2.1 meter cord
- 1 – 9 volt, medical grade, international power supply with 6 foot / 1.8 meter cord and universal plugs
- 1 – Handpiece surface mounting bracket with double stick adhesive tape
- 1 – VALO Barrier Sleeves sample pack
- 1 – Amber colored safety glasses

Overview of Controls:

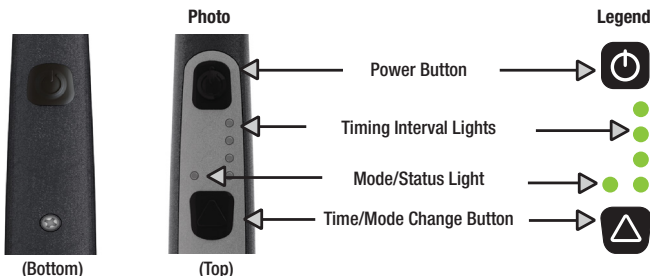Instructions for Use:

- 1 - Connect the 9-volt power cord to the handpiece cord.
- 2 - Plug the power cord into any electrical outlet (100 -240 VAC). The timing lights will illuminate indicating the light is ready for use. The VALO GRAND Corded curing light handpiece beeps twice when powering on.
- 3 - Place the VALO GRAND Corded curing light handpiece into a standard dental unit mounting bracket or accessory mounting bracket until ready to use.
- 4 - Prior to use, place barrier sleeve over VALO GRAND Corded curing light wand.
- 5 - Select the desired mode: VALO GRAND Corded curing light has 3 modes: Standard Power, High Power Plus, and Xtra Power. Each mode is identified by the Mode/Status Light (green = Standard Power, orange = High Power Plus, flashing orange = Xtra Power). To change modes, hold the Time/Mode Change button for 2 seconds and release. The Mode/Status light will change to the next mode.

NOTE: VALO GRAND Corded curing light is programmed to cycle from the Standard Power to the High Power Plus to the Xtra Power mode in sequence. For example, to change from the Standard Power mode to the Xtra Power mode, it is necessary to cycle into the High Power Plus mode and then to the Xtra Power mode.

VALO GRAND Corded curing light always stores its last used timing interval in each mode and will default back to that timing interval whenever the modes are changed, or if the unit is unplugged and loses power.

SLEEP MODE: VALO GRAND Corded curing light will go into SLEEP mode after 1 hour of inactivity, as indicated by the green mode light blinking slowly. Press any button to wake VALO GRAND Corded curing light and begin curing. VALO GRAND Corded curing light will automatically return to the last setting used.

Precaution: Always wear amber colored safety glasses while the curing light is illuminated.

Quick Mode Guide


Mode	Standard Power				High Power Plus				Xtra Power
Power Button									
Mode/Timing LEDs									
Time Button									
Time Options	5s	10s	15s	20s	1s	2s	3s	4s	3s only
To Change Time	Press and release Time Button quickly to cycle through time options.								
To Change Modes	Press and hold the Time Button for 2 seconds and release. VALO GRAND Corded curing light will cycle to next Mode.								
Legend	Solid LEDs Blinking LEDs								

Quick Curing Guide: Recommended Curing Times for Optimal Results with VALO GRAND Corded curing light

Exposure times may need to be adjusted due to composite reactivity, shade, distance from the light lens to the composite, and depth of composite layer if it is over 2mm.

Mode	Standard Mode	High Power Plus Mode	Xtra Power Mode
Per Layer	One 10 second cure	Two 4 second cures	One 3 second cure
Final Cure	Two 10 second cures	Three 4 second cures	Two 3 second cures

Quick Warning Guide

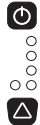Warnings	
Call Customer Service for Repair	Call Customer Service for Repair
No sound. Flashing, 2 seconds Allows operation	Continuous 3 beeps Prohibits operations
	

CURING MODE: Standard Power mode

USES: Curing of restorative materials with photo initiators.

TIMING INTERVALS: 5, 10, 15, 20 seconds.

- VALO GRAND Corded curing light defaults to this mode when it is INITIALLY powered on. The Status Light is green and the green Timing Lights are solidly illuminated.
- To change timing intervals, quickly press the Time/Mode Button.
- Press the Power Button to cure. To stop curing prior to completion of a timing interval, press the Power Button again.

CURING MODE: High Power Plus mode

USES: Curing of restorative materials with photo initiators.

TIMING INTERVALS: 1, 2, 3, 4 seconds.

- From Standard Power mode, press and hold the Time/Mode Change Button for 2 seconds. The green Timing Lights will illuminate and flash. The Status Light will illuminate as a solid orange light, indicating High Power Plus mode. A two-second power tack is the most commonly used timing interval in this mode.
- To change timing intervals, quickly press the Time/Mode Button.
- Press the Power Button to cure. To stop curing prior to the completion of a timing interval, press the Power Button again.
- To return to Standard Power mode, press and hold the Time/Mode Change Button for 2 seconds, this will cycle to Xtra Power Mode, then release, hold for 2 seconds, and release. The green light will be solid with solid green timing lights, indicating Standard mode.

CURING MODE: Xtra Power mode

USES: Curing of restorative materials with photo initiators.

TIMING INTERVAL: 3 seconds only (Note: there is a 2 second safety delay at the end of each curing cycle).

- From the Standard Power mode, press the Time/Mode Change Button for 2 seconds, release, press again for 2 seconds, and release. Three of the green Timing Lights and the orange Status Light will illuminate and flash, indicating Xtra Power mode.
- Press the Power Button to cure. To stop curing prior to completion of a timing interval, press the Power Button again.
- To return to the Standard Power mode, press and hold the Time Button for 2 seconds.

WARNING: In Xtra Power mode, DO NOT expose soft tissue for more than 10 seconds. The Xtra Power mode has a 2-second safety delay to limit heating during consecutive curing. At the end of the delay, beeping indicates unit is ready for continued used. If a longer cure is needed, allow 10 seconds between consecutive cures or consider a dual-cure product (composite or adhesive).

Note: If it is necessary to cure near soft tissue, cooling with an air stream is an option for keeping the area cooler.

Duty Cycle: The VALO GRAND Corded curing light is designed for short-term operation. At maximum ambient temperature (32°C) 1 minute ON back-to-back cycling, 30 minutes OFF (cooling-off period).

Accessory Bracket Instructions

- 1- Bracket should be mounted to a flat, oil-free surface.
- 2- Clean surface with rubbing alcohol.
- 3- Peel backing off of the bracket's adhesive tape.
- 4- Position bracket so the handpiece lifts upward when removed. Press firmly into place.

Maintenance and Cleaning

VALO GRAND Corded curing light is a sealed unit with a hard surface and a scratch resistant glass lens. After each use, moisten a gauze or soft cloth with an approved surface disinfectant and wipe the surface and lens.

Routinely check the lens for cured dental resins. If necessary, use a non-diamond dental instrument to carefully remove any adhered resin.

The barrier sleeve helps prevent cross contamination and helps keep dental composite material from adhering to the surface of the lens and wand body.

Light meters differ greatly and are designed for specific light guide tips and lenses. Ultradent recommends checking VALO GRAND Corded curing light in Standard Power mode. NOTE: the true numeric output will be skewed due the inaccuracy of common light meters and the custom LED pack VALO GRAND Corded curing light uses.

PRODUCT SAFETY WARNINGS:

- DO NOT autoclave.
- DO NOT immerse in any kind of ultrasonic bath, disinfectant, cleaning solutions, or any kind of liquid.
- PRECAUTION: This unit may be susceptible to strong magnetic or static electric fields, which could disrupt the programming. If you suspect this has occurred, unplug the unit momentarily and then replug it into the outlet.
- DO NOT wipe down the VALO GRAND Corded curing light with caustic or abrasive cleaners. See lists of acceptable cleaners below:

ACCEPTABLE CLEANERS:

- Lysol Brand III Disinfectant Spray (Recommended)
- Isopropyl alcohol
- Ethyl alcohol-based cleaners
- Lysol®* Concentrate (alcohol-based only)
- Cavicide™** products (Non-Bleach)**

DO NOT USE:

- Strong alkali detergent of any type, including hand soaps and dish soaps
- Bleach-based cleaners (e.g. - Clorox™*, Sterilox™**)
- Hydrogen Peroxide based cleaners
- Abrasive Cleansers (e.g. - Comet Cleanser™**)
- Acetone or hydrocarbon based cleaners
- MEK (Methyl Ethyl Ketone)
- Birx®*
- Gluteraldehyde
- Quaternary Ammonium Chloride salt-based cleaners (except Cavicide)
- Cavicide1 solution or wipes

* Trademark of a company other than Ultradent.

** If used, it may fade the color


Troubleshooting Guide

If the solutions suggested below do not rectify the problem, please call Ultradent at 800.552.5512.

Outside the United States, call your Ultradent distributor or dental dealer.

Repairs are only to be performed by authorized service personnel. Ultradent to provide service personnel with documentation to perform repairs.

Problem	Possible Solutions
Light will not turn on	1. Press the Time/Mode or either Power Button to wake from Power Save Mode. 2. Check that both cords are firmly connected together and to the electrical outlet. 3. Confirm power to the wall outlet.
Light does not stay on for desired time	1. Check Mode and Timing lights for correct time input. 2. Confirm all cord connections are fully seated. 3. Unplug and re-plug power cord into the electrical receptacle.
Light is not curing resins properly	1. Check lens for residual cured resins/composites. (see "Maintenance and Cleaning") 2. Using proper amber UV eye protection, verify the LED lights are working. 3. Check power level with light meter. (See "Maintenance and Cleaning.") 4. Check expiration date on curing resin. 5. Ensure proper technique is being followed (adhesive/composite) by manufacturer recommendations.

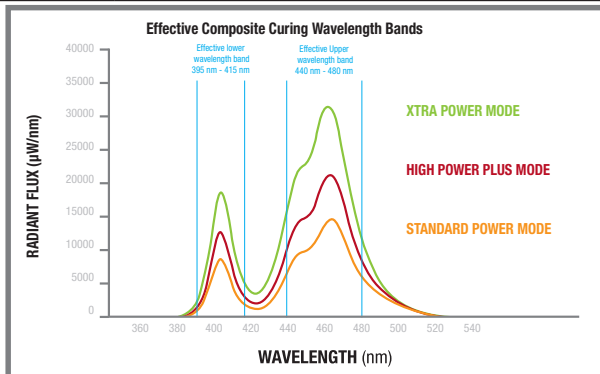Technical Information	VALO GRAND Corded LED Curing Light			
Lens	Diameter 11.7 mm			
Wavelength range	385nm – 515nm (see qualification below) Effective output Power of VALO GRAND Corded curing light falls within the following wavelength range: • 395nm <= EP <= 480nm Minimal and insignificant power can be found in wavelength ranges: • 385nm – 395nm and 480nm – 515nm			
Light Intensity Table	Nominal Radiant Exitance Comparison Chart			
	Measurement Instrument	* Demetron L.E.D. Radiometer	* MARC spectrum analyzer	‡ Gigahertz spectrum analyzer
	Aperture of Meter	7 mm	3.9 mm	15 mm 15 mm
	Standard Power (±10%)	1000 mW/cm²	1200 mW/cm²	900 mW/cm² 970 mW
	High Power Plus (±10%)	1600 mW/cm²	1800 mW/cm²	1500 mW/cm² 1615 mW
	Xtra Power (±20%)	N/A	3200 mW/cm²	2100 mW/cm² 2260 mW

Radiant Exitance will vary based on instrument capability, measurement method and light placement.

* Demetron radiometers and MARC spectrum analyzers should be used as reference only due to having smaller apertures than VALO GRAND Corded curing lights.

* Demetron radiometers should be used as reference only due to limitations in power and spectral response.

‡ Radiant Exitance conforms to ISO 10650 when measured with a Gigahertz spectrum analyzer.

Wand	Weight - 8 ounces/226 grams (with cord) Length - 9.26 inches/23.5 cm Breadth -.74 inches/1.9 cm	Width - .79 inches/ 2 cm Cord length - 7 feet/2.1 meters Ratings - CE, IEC 60601
Power Supply	Output - 9VDC at 2A Input - 100VAC to 240VAC Ultradent P/N 5930 VALO GRAND Corded curing light Power Supply with Universal Plugs	Rating - Medical Grade Class II (UL, CE, RoHS, WEEE) Cord length - 6 feet/1.8 meters VALO GRAND Corded curing light power supply provides isolation from MAINS power
Operating Conditions	Temperature: +10°C to +32°C (+50°F to +90°F) Relative Humidity: 10% to 95% Ambient Pressure: 700 hPa to 1060 hPa	
Storage and Transport Conditions:	Device Storage and Transport: • Temperature: +10°C to +40°C (+50°F to +104°F) • Relative Humidity: 10% to 95% • Ambient Pressure: 500 hPa to 1060 hPa	

Guidance and Manufacture's Declaration for Electromagnetic Emissions

The VALO GRAND Corded curing light is intended for use in the electromagnetic environment specified below. The customer or user should ensure that it is used in such an environment.

WARNING: Use only authorized accessories, cables, and power supplies to prevent improper operation, increased electromagnetic emissions or decreased electromagnetic immunity.

Emissions Test	Compliance	Electromagnetic environment - guidance
RF emissions: CISPR 11	Group 1	VALO GRAND Corded curing light uses a Globtek medical grade 9VDC adaptor, operates with brown-out protection, and provides limited EMI, RF, and surge suppression. VALO GRAND Corded curing light uses electrical and electromagnetic energy only for their internal functions. Therefore, any RF emissions are very low and are not likely to cause interference in nearby electronic equipment. The VALO GRAND Corded curing light is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings for domestic use.
RF emission: CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	COMPLIES	

Guidance and Manufacture's Declaration for Electromagnetic Immunity

The VALO GRAND Corded curing light is intended for use in the electromagnetic environment specified below. The customer or user should ensure that it is used in such an environment.

IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 15 kV air	± 8 kV contact ± 15 kV air	Physical environment should be restricted to the following: 1. IP Code: IP20 2. Do not immerse in liquid. 3. Do not use around flammable gas. Unit is Non-APG and Non-AP. 4. Storage humidity range: 5% - 95% 5. Storage temperatures range: 10° C - 40° C Mains power quality should be that of a typical residential, commercial, or hospital, or military environment
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines Note 1: VALO GRAND Corded curing light has no I/O ports	
Surge IEC 61000-4-5	± 1 kV line to line ± 2 kV line to earth	± 1 kV line to line ± 2 kV line to earth	

Voltage, dips, shorts, interruptions and variations on the power supply input lines IEC 61000-4-11	<5% U (>95% dip in U for 0.5 cycle) 40% U (60% dip in U for 5 cycles) 70% U (30% dip in U for 25 cycles) <5% U (>95% dip in U for 5 s)	<5% U (>95% dip in U for 0.5 cycle) 40% U (60% dip in U for 5 cycles) 70% U (30% dip in U for 25 cycles) <5% U (>95% dip in U for 5 s) Note 2: Self recovers	Mains power quality should be that of a typical residential, commercial, hospital, or military environment. The Globtek 9VDC medical grade adaptor that is supplied with the VALO GRAND Corded curing light operates from mains ranging from 100VAC – 240VAC and is capable of limited brown out, EMI, and surge protection. If the VALO GRAND Corded curing light user requires continued operations without mains interruption, or the mains in any particular region of a country are considered bad due to continuous brown-out, black-out, or excessively noisy power conditions, it is recommended that the VALO GRAND Corded curing light be powered from an uninterruptible power supply or the customer purchase a VALO GRAND Corded curing light unit.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical, residential, home health care, commercial, hospital, or military environment.

NOTE: U is the A.C. mains voltage prior to application of the test level

Note 1: The VALO GRAND Corded curing light is not equipped with any ports or any accessible I/O lines.

Note 2: If there is a drop of 95% in Mains voltage VALO GRAND Corded curing light will not operate. It has no internal energy storage mechanism. The VALO GRAND Corded curing light will turn off. When power levels are restored, VALO GRAND Corded curing light will restart and return to the same state before power loss. The VALO GRAND Corded curing light will self-recover.

Guidance and Manufacture's Declaration for Electromagnetic Immunity for non-life support systems

The VALO GRAND Corded curing light is intended for use in the electromagnetic environment specified below. The customer or user should ensure that it is used in such an environment.

IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Conduction RF	3 Vrms	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the VALO GRAND Corded curing light, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = \left[\frac{3,5}{V_1} \right] \sqrt{P}$ $d = \left[\frac{3,5}{E_1} \right] \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = \left[\frac{7}{E_1} \right] \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,* should be less than the compliance level in each frequency range.^b Interference may occur in the vicinity of equipment marked with the following symbol: 
IEC 61000-4-6	150 kHz to 80 MHz	150 kHz to 80 MHz	
Radiated RF	3 V/m	3 V/m	
IEC 61000-4-3	80 MHz to 2.5 GHz	80 MHz to 2.5 GHz	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the VALO GRAND Corded curing light is used exceeds the applicable RF compliance level above, the VALO GRAND Corded curing light should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the VALO GRAND Corded curing light.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Guidance and Manufacturer's Declaration for recommended separation distances between portable and mobile RF communications equipment and the VALO GRAND Corded curing light

WARNING: Portable RF communications equipment may degrade performance if used closer than 30 cm (12 in").

The VALO GRAND Corded curing light is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The user of the VALO GRAND Corded curing light can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the VALO GRAND Corded curing light as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (P in Watts)	Separation distance according to frequency of transmitter (meters)		
	150 kHz – 80 MHz	80 MHz – 800 MHz	800 MHz – 2.5 GHz
	$d = \left[\frac{3,5}{V_1} \right] \sqrt{P}$	$d = \left[\frac{3,5}{E_1} \right] \sqrt{P}$	$d = \left[\frac{7}{E_1} \right] \sqrt{P}$
0.01	0.12 meters	0.035 meters	0.07 meters
0.1	0.37 meters	0.11 meters	0.22 meters
1	1.7 meter	0.35 meter	0.7 meter
10	3.7 meter	1.11 meter	2.22 meter
100	11.7 meters	3.5 meters	7.0 meters

The VALO GRAND Corded curing light has been tested according to IEC 60601-1-2:2014 and passed under radiated field strengths of 10 V/m between 80-MHz to 2.5 GHz. The value of 3Vrms corresponds to V1 and the value 10 V/m corresponds to E1 in the formulas above.

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Warranty for VALO GRAND Corded curing light units:

Ultradent hereby warrants that this instrument shall, for a period of 3 years**, conform in all material respects to the specifications therefore as set forth in Ultradent's documentation accompanying the product and be free from any defects in materials/or workmanship. This warranty applies solely to the original purchaser and is not transferable. All defective products are to be returned to Ultradent. There are no user service components of the VALO GRAND Corded curing light system. Tampering with VALO GRAND Corded curing light will void its warranty.

The VALO GRAND Corded warranty does not cover customer caused damage. For example; if a VALO GRAND Corded is misused or dropped and the lens breaks, the customer would be responsible to pay for any necessary repairs.

**With sales receipt indicating the date of sale to the dentist.

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© 2019 Ultradent Products, Inc. All Rights Reserved.
Manufactured by Ultradent Products, Inc. 505 West Ultradent Drive
Ultradent Products GmbH