




SOOOO SUPER
WITH **FUTAR[®]**

BITE REGISTRATION

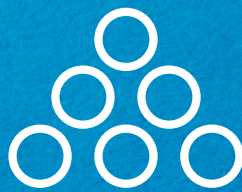

KETTENBACHDENTAL

Simply intelligent

FUTAR® – SIMPLY IMPRESSIVE!

Bite registration

The **Futar®** product family includes six different bite registration materials based on A-silicone for making accurate occlusal records.

SUPER IN QUALITY AND PROGRESS

- "Made in Germany" with its own research and development as well as production in Germany

SUPER THROUGH VERSATILITY

- Choose between high and highest final hardness in addition to different processing times!

SUPER FOR OVER 20 YEARS

- Trust Germany's market leader*!

Do not make compromises with your bite registration!



FUTAR® – THE ADVANTAGES AT A GLANCE

Futar® Cut & Trim Fast:

- Economical thanks to the **use of the shorter yellow mixing tips** (up to three more registrations per cartridge)
- Extra-high final hardness (Shore-D 35)
- Ideal for processing with a mill and / or scalpel
- Scannable for using with CAD/CAM in your practice!



Futar® D / Futar® D Fast:

- Extra-high final hardness (Shore-D 43)
- Offers very good millability
- Also available with fast-setting characteristics: **for more patient comfort!**



The Futar® family can
do even more:

Futar® / Futar® Fast:

- High final hardness (Shore-A 90)
- Hard in setting, flexible in processing: easy to cut
- Also available with fast-setting characteristics: **for more patient comfort!**



Futar® D Slow:

- Extra-high final hardness (Shore-D 43)
- Extra-long processing time for time-consuming applications, e.g. for myocentric bite registration, fixation of impression posts or individual contouring of the margins
- Versatile use also in the laboratory: as insulating material or for the fabrication of model segments








- Accurate results thanks to excellent flow properties
- Low elasticity reduces springiness in the model assignment
- Surplus can be broken off easily

**SIMPLE,
CLEAN AND
FAST TO USE!**

FUTAR® – PERFORMANCE RUNS IN THE FAMILY








Dental arch

 Bite registration material	 Mixing	Working time at 23 °C / 74 °F ≤	Intraoral setting time at 35 °C / 95 °F ≥	Total setting time* ≥	Hardness (approx.) Shore
Futar®		30 sec.	90 sec.	2:00 min.	A 90
Futar® D		30 sec.	90 sec.	2:00 min.	D 43
Futar® D Slow		90 sec.	3:00 min.	4:30 min.	D 43

*Total setting time (removal time from mouth) from the start of mixing

Segment

 Bite registration material	 Mixing	Working time at 23 °C	Intraoral setting time	Total setting time* ≥	Hardness (approx.) Shore
Futar® Fast		15 sec.	45 sec.	1:00 min.	A 90
Futar® D Fast		15 sec.	45 sec.	1:00 min.	D 43
Futar® Cut & Trim Fast		15 sec.	45 sec.	1:00 min.	D 35

*Total setting time (removal time from mouth) from the start of mixing

**CREATE
YOUR BITE
REGISTRATION
IN JUST ONE
MINUTE!**

FUTAR® – THE RANGE

Futar®	Futar® Fast	Futar® D	Futar® D Fast
Intro pack	50 mL cartridge, 6 green mixing tips		
Normal pack	2 × 50 mL cartridges, 6 green mixing tips		
Bonus pack	10 × 50 mL cartridges, 30 green mixing tips		

Futar® Cut & Trim Fast	
Intro pack	50 mL cartridges, 8 yellow mixing tips
Normal pack	2 x 50 mL cartridges, 8 yellow mixing tips
Bonus pack	10 x 50 mL cartridges, 40 yellow mixing tips

Futar® D Slow	
Normal pack	2 × 50 mL cartridges, 6 green mixing tips



Kettenbach GmbH & Co. KG
Im Heerfeld 7
35713 Eschenburg, Germany
Phone +49 (0)2774 7050
info@kettenbach.com
www.kettenbach-dental.com


KETTENBACHDENTAL
Simply intelligent

022744_0919

